

# Starkey's Complete Style Portfolio



Genesis AI marks a new beginning for the hearing industry.

## GENESIS<sup>AI</sup>



2.4 GHz HEARING AIDS

Feature/ Product Style	CIC 312	ITC/HS R	ITER	mRIC R	RIC RT	RIC 312
Battery Type	312	Rechargeable Lithium-ion	Rechargeable Lithium-ion	Rechargeable Lithium-ion	Rechargeable Lithium-ion	312
Telecoil			•		•	•
CROS System					•	•
Tinnitus Technology	•	•	•	•	•	•

\*IIC NW excluded from 2.4 GHz platform

2.4 GHz ACCESSORIES  
(StarLink)



### Remote Control 2.0

Redesigned for maximum simplicity, this remote lets users easily control volume, change programs and quickly access a favorite feature. Comes with replaceable cell battery and handy lanyard.



### TV Streamer

Offering excellent sound quality, stream audio from the television or other electronic audio source directly to hearing aids.



### Remote Mic Plus

Use as a remote microphone, or easily stream binaurally from four different audio sources: Bluetooth®, Telecoil, FM, or Line-In.



### Mini Remote Mic

Enjoy one-on-one conversations in noisy environments with our small, easy-to-use Mini Remote Microphone.



### Table Microphone

With 8 built-in microphones, users can enjoy group settings like family gatherings or noisy restaurants.

CHARGERS



### StarLink Premium Charger 2.0 (RIC RT)



### StarLink Premium Charger 2.0 (mRIC)



### StarLink Premium Charger 2.0 (Custom)


Our StarLink Premium Charger 2.0 is the ROES default for all Genesis AI rechargeable hearing aid orders. This charger includes an on-board lithium-ion battery that provides up to three full charges on-the-go. The tailored interior charges RICs and custom hearing aids easily and securely. It has internal LED charging lights for the charger battery, as well as each hearing aid battery that show charging in progress and charging complete. The Premium charger includes on-board turbo charging for three hours of hearing aid use time in just 10 minutes of charging.

*StarLink 2.0 Non-premium chargers are available to order in ROES as a selectable option with new rechargeable orders. Non-premium chargers do not include an on-board battery and must be plugged into an external power source (i.e., wall outlet) to function/charge the hearing aids.*

*Standalone StarLink Premium and Non-premium chargers are available to order through Parts and Accessories.*



# My Starkey App

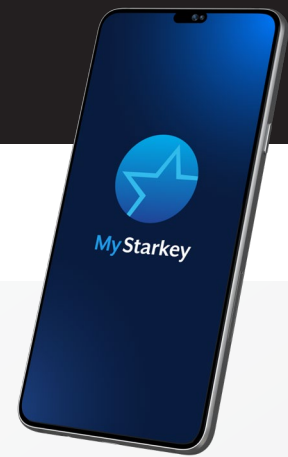

## Personal Health Monitoring Technology Features:

Available with **Genesis AI** products



### Edge Mode+

Edge Mode+ is now automatic. It continues to scan, adapt, and optimize as the acoustic environment changes.



### Hearing Health and Activity Tracking

Monitors vigorous physical activity and personal health monitoring features like active listening.



### Patient Education

A complete library of fresh "how-to" content and videos.



### Smart Assistant

Voice-activated audio adjustments, reminders, answers to hearing aid topics, and more.



### Hearing Control

Change volume and switch memories with your smartphone.



### Personalized Programs

Create custom programs and geotag programs so adjustments occur automatically when you return to that location.



### Self Check

Allows you to perform a diagnostic test of your hearing aid components.



### Fall Alert

With integrated sensors, Fall Alert can detect falls and send alert messages to selected contacts.



### Social Interaction Updates

Can now distinguish patient's voice vs. others' to more accurately determine time spent in conversation.



### Personalized Control

Make personalized, real-time adjustments.



### 2-Way Audio

Take and end calls directly from your hearing aids instead of your Apple® device. Plus, hearing aid microphones pick up and stream your voice directly back to iPhone and iPad\*, allowing you to talk hands-free.



### Language Translation

Communicate with people who speak other languages via the easy-to-use translation tool in your My Starkey app (requires an internet connection).



### Tap Control

Tap your ear twice to start & stop streaming audio from a 2.4 GHz accessory.



### TeleHear

Remote Programming allows for virtual adjustments to hearing aids.



### Find My Phone

Using hearing aids and Smart Assistant, you can now locate your lost or misplaced phone by having it ring. No other hearing aid offers this capability.



### Mask Mode

A custom memory that boosts speech audibility.



### Auto Experience Manager

Automatically adjusts loudness over time to help you transition your new hearing experience comfortably.



### Intelligent Reminders

Hearing Aid Maintenance reminds patients to run diagnostic self-check on their hearing aids.



### Phone Calls

Answer calls and hear the conversation streamed directly to your hearing aids using compatible devices.



### Transcribe

Conversations are transcribed so you can read what's being said into text that can be saved or shared.



### Voice AI\*

This always-learning speech enhancement feature uses an iOS smartphone to deliver an unprecedented experience for those with moderate loss of 50 dB HL or greater.



### Find My Hearing Aids

Easily locate lost hearing aids.



### Adaptive Car Mode

Hearing aids automatically change settings to reduce the annoying sounds of driving.



### Reminders

Personalized or on-screen reminders for daily tasks such as taking medications.



### Siri Integration

Have Siri® read texts and emails directly to your hearing aids (available with Apple devices).



### Adjust Sound

When streaming "Adjust Sound" appears on app Home screen for access to EQ controls.

\*Requires an iOS device.

The app may look slightly different depending on your phone.

Starkey, the Starkey logo, Genesis, StarLink and Edge Mode are registered trademarks of Starkey Laboratories, Inc. Google Play is a trademark of Google, LLC. | The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks is under license. | Siri, iPhone, iPad and Apple are registered trademarks of Apple Inc.

©2024 Starkey Laboratories, Inc. All Rights Reserved. P00001465 4/24 FLYR4110-02-EN-GS

